

vygon**vet**  
PRODUCT FOCUS

MY VYGON VET PET WINNER 2023



# WHY CHOOSE VYGON?




OVER  
**60**  
YEARS  
EXPERIENCE

**PATIENT CARE**  
LIES AT THE HEART  
OF WHAT WE DO



**CLINICAL**  
EXPERTISE

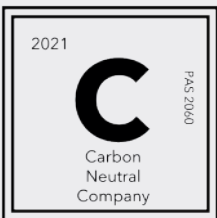

WE ARE A  
**CARBON**  
**NEUTRAL**  
COMPANY

MORE THAN  
**100 MILLION**  
PRODUCTS IN  
**110 COUNTRIES**



# CONTENTS

Why choose Vygon? \_\_\_\_\_ 2

About Vygon Vet \_\_\_\_\_ 3

Our commitment to the environment \_\_\_\_\_ 3

Needle-free devices \_\_\_\_\_ 4

IV management \_\_\_\_\_ 8

CVCs \_\_\_\_\_ 10

Nerve blocks \_\_\_\_\_ 12

Cannula \_\_\_\_\_ 14

Mostcare Up \_\_\_\_\_ 15

VYSET procedure packs \_\_\_\_\_ 16

Oncology \_\_\_\_\_ 18

Enteral feeding \_\_\_\_\_ 20

Theatre \_\_\_\_\_ 22

Urinary \_\_\_\_\_ 25

Aitecs infusion pumps \_\_\_\_\_ 26

# OUR SERVICE OFFERING



## Customer Service & Technical Support

Talk through your enquiry with our  
dedicated Teams.



## Education & Training

We offer our customers a variety of  
valuable and comprehensive training  
options to help you and your teams  
meet your training requirements.

This Product Focus features a  
selection of our products.  
Visit [vygonvet.co.uk](http://vygonvet.co.uk) for our full range.



A leading supplier of single-use veterinary consumables for more than 20 years, Vygon Vet has built an enviable reputation for high quality products and excellent customer service.

Working closely with practice teams focused on delivering the best care for their patients, our priority is to tailor our products and support to ensure surgical procedures run smoothly and efficiently and treatments and therapies deliver faster and better recovery. This philosophy is behind the success of all our products including the popular needle-free IV access system Bionector, atraumatic long term feline urinary catheter Katkath and Vyset custom packs.

In addition, we also focus on supporting practices to push through treatment boundaries and benefit from innovation. Ground-breaking products include the Plexygon Nerve Stimulator, to access the patient benefits of regional anaesthesia, and Qimono, our closed system transfer device, which provides protection for veterinary teams whilst delivering life-saving chemotherapy.

Vygon Vet is a division of leading healthcare and medical products supplier Vygon (UK) Ltd.

## OUR COMMITMENT TO THE ENVIRONMENT

Vygon UK has been working towards creating a sustainable future for many years and we are a Carbon Neutral Company. We believe that sustainability isn't just about meeting current needs, but more importantly, it is about ensuring we are here for the long term and are paving the way for a bright tomorrow. We strive to deliver enhanced Corporate Social Responsibility (CSR), ensuring we manage the social, economic, and environmental effects of Vygon UK's operations responsibly in line with public expectations.




DISCOVER OUR JOURNEY  
scan with your smart device



# Needle-free devices

## Single bionector

Code	Description	Box qty
0896.03	bionector closed system needle-free device	50

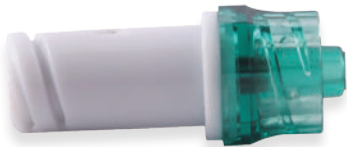

## bionector T-Piece

Allows intermittent injection via latex-free membrane during simultaneous infusion.

Code	Description	Box qty
1VTUSBN-448	10cm luer-lock T-piece extension connector with bionector	50



## bionector T-piece with two needle-free access ports

Totally needle-free T-piece device, reduces risk of needle-stick injuries to both patients and staff, as well as reducing the risk of infection.

Code	Box qty
1VTUSBN-450	50



## bionector vial caps

For use when aspirating drugs from capped vial and drug reconstitution.

Code	Size (mm)	Box qty
0821.01	14	50
0821.02	20	50




## bionector and bag spike

Code	Box Qty
0856.11	50





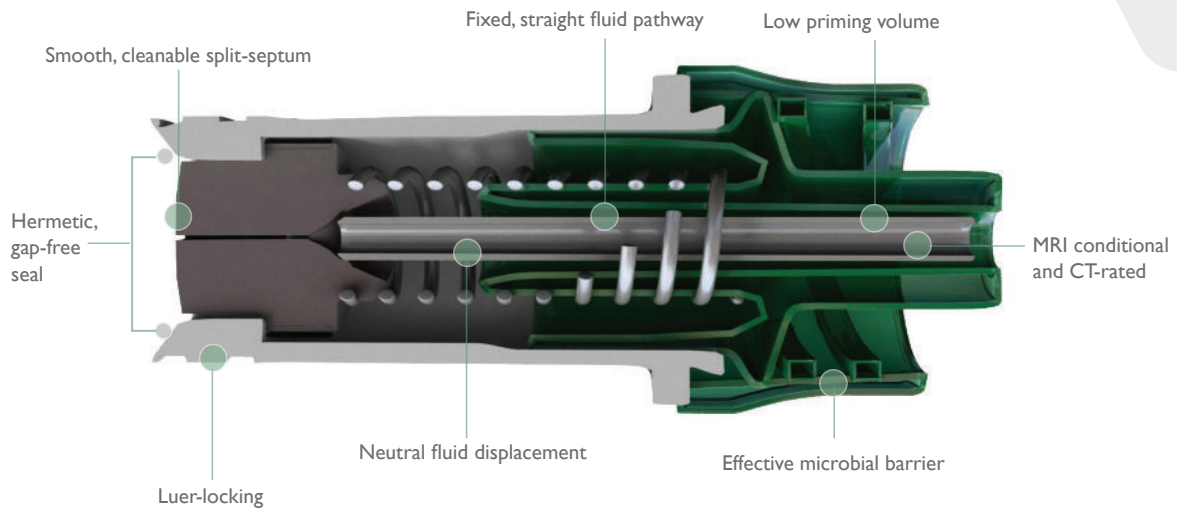

## Exceeding global opinion leaders' recommendations

An established market-leading needle-free device which meets the full range of global opinion leaders' recommendations for reducing CRBSIs. It has been proven to provide an effective barrier against microbial ingress and help standardise practice by combining a fixed, straight fluid pathway with innovative neutral displacement technology.

## MRI conditional and CT-rated

Bionector is proven not to represent any risk to either patients or practitioners during an MRI of up to three Teslas. CT-rated for use with power injectors, Bionector has a maximum pressure resistance of 350psi and a maximum flow rate of 10mL/s.

## Neutral displacement

Leading the way with a neutral fluid displacement, a specific post-flushing clamping sequence is not required. This helps prevent blood reflux and reduce catheter occlusions.

## Clinically proven

Backed up by a robust library of clinical studies Bionector is proven to be easy to clean and clear. Its smooth split-septum fits tightly into the device housing ensuring it is free from any gaps. The straight, fixed fluid pathway has been proven 'easy to clear', designed to provide the most direct and least tortuous route with no moving parts (such as mechanical valves), which reduces the surface area available for biofilm formation.

## Low deadspace

The straight fluid pathway is proven 'flushable' for macro and microscopic particles such as blood. This is due to a minimal deadspace of just 0.018mL allowing for a low flushing volume (5mL) to clear the device.



# Needle-free devices

## bionector octopus double lumen

Octopus needle-free extension set, including double lumen Bionector octopus extension with 10cm tube and clamp PUR.

Code	Length (cm)	Flow rate (mL/min)	ID (mm)	OD (mm)	Box qty
0841.264	10	100	1.5	2.5	50



## bionector octopus double lumen with two anti-reflux valves

Octopus needle-free extension set, including double lumen Bionector octopus extension and two anti-reflux valves with 0.4mL capacity.

Code	Length (cm)	Flow rate (mL/min)	ID (mm)	OD (mm)	Box qty
0841.232	3	120	1.5	2.5	10



## bionector octopus triple lumen

Octopus needle-free extension set, including triple lumen Bionector octopus extension with 10cm tube and clamp PUR.

Code	Length (cm)	Flow rate (mL/min)	ID (mm)	OD (mm)	Box qty
0841.364	10	100	1.5	2.5	10



## bionector octopus triple lumen with three anti-reflux valves

Octopus needle-free extension set, including triple lumen Bionector octopus extension and three anti-reflux valves with 0.4mL capacity.

Code	Length (cm)	Flow rate (mL/min)	ID (mm)	OD (mm)	Box qty
0842.312	3	120	1.5	2.5	10




## 3M™ Curost™ disinfecting port protectors

A disposable passive disinfection device. Curost port protectors ensures needle-free access devices are always clean, therefore removing the need to actively disinfect before and after each access.

Code	Description	Box qty
VWJIM001-270	Curost™ disinfecting port protector compact patient box	30

Curost™ is a trademark of 3M Company.

**3M**  
Authorised Distributor  
Medical Products

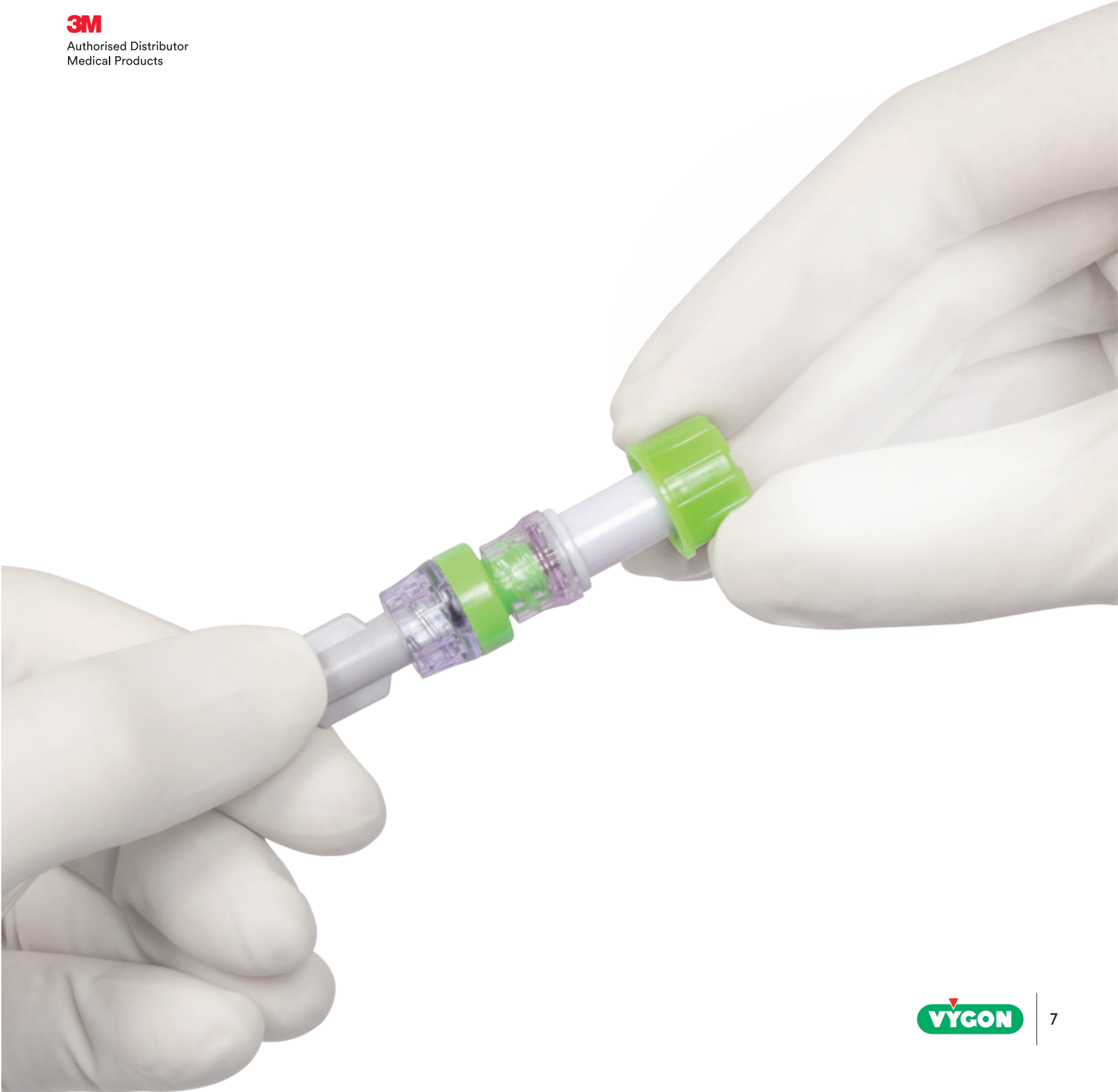


# IV management

## lectrocath

lectrocath 50cm micro bore male/female luer-lock PE extension set.

Code	Length (cm)	Capacity (mL)	ID (mm)	OD (mm)	Box qty
1155.05	0.7	2.0	50	1.0	50

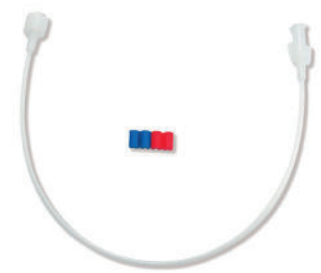

## lectrospiral








Coiled IV extension line in polyurethane with luer-lock connectors.

Code	Length (cm)	Capacity (mL)	ID (mm)	OD (mm)	Box qty
1155.60	100	1.2	1.0	2.0	15
1155.65	150	1.6	1.0	2.0	15
1155.70	200	2.0	1.0	2.0	15
1155.80	300	2.9	1.0	2.0	15
1155.90	400	3.8	1.0	2.0	15
1159.60	100	2.0	1.5	2.5	15
1159.65	150	2.8	1.5	2.5	15
1159.70	200	3.6	1.5	2.5	15
1159.80	300	5.2	1.5	2.5	15
1159.90	400	6.9	1.5	2.5	15



## microflex

Silicone-coated needle with double wings.

Code	Colour	Size (G)	OD (mm)	Length (cm)	Box qty
000246044		27	0.4	30	100
000246054		25	0.5	30	100
000246064		23	0.6	30	100
000246074		22	0.7	30	100
000246084		21	0.8	30	100
000246094		20	0.9	30	100
000246114		19	1.1	30	100






v-green IV extension line

PE extension line with male to female luer-lock connections. Narrow bore to minimise priming volume, ideal for drug infusion in animals of all sizes.

Code	Length (cm)	Capacity (mL)	ID (mm)	OD (mm)	Box qty
7.1100.05	50	0.7	1.0	2.0	50
7.1100.10	100	1.2	1.0	2.0	50
7.1100.15	150	1.6	1.0	2.0	50
7.1100.20	200	2.0	1.0	2.0	50
7.1101.05	50	1.1	1.5	2.5	50
7.1101.10	100	2.0	1.5	2.5	50
7.1101.15	150	2.8	1.5	2.5	50
7.1101.20	200	3.6	1.5	2.5	50

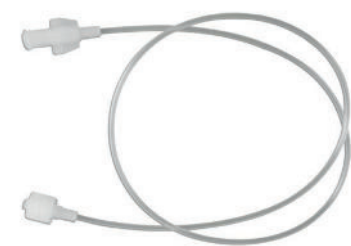

vyclic three-way stopcock

Lipid resistant, lockable three-way stopcock in PC-PE, with one rotating male luer and two female luer-lock ports. Colour-coded plugs for venous or arterial identification, and the handle position indicates open ports.

Code	Description	Box qty
087620	Lockable three-way stopcock in PC-PE	50



protect-a-line extension line

Code	Length (cm)	Capacity (mL)	ID (mm)	OD (mm)	Box qty
0832.01	150	1.4	1.0	2.5	50
0832.211	200	1.3	1.5	3.0	50




# CVCs

## multicath forward

Radiopaque multi-lumen central venous catheter with flexible tip. 1cm graduations start at 9cm from the distal tip. Each lumen has a different hub colour to aid identification. The guidewire feeder aids insertion whilst still allowing the user to sense any resistance. Each kit includes an introducer needle and a short PTFE cannula to give the user clinical choice for insertion.


Contents:

- 1 x multicath forward catheter
- 1 x Introducer needle 18G x 70mm
- 1 x Mosquito introducer cannula 18G x 64mm
- 1 x "J" guidewire in feeder; length 60cm, OD 0.9mm
- 1 x Vessel dilator
- 2 x Latex-free injection membranes
- 1 x Mini scalpel N
- 1 x 5mL luer-slip syringe
- 1 x Secondary fixation wing and anchor

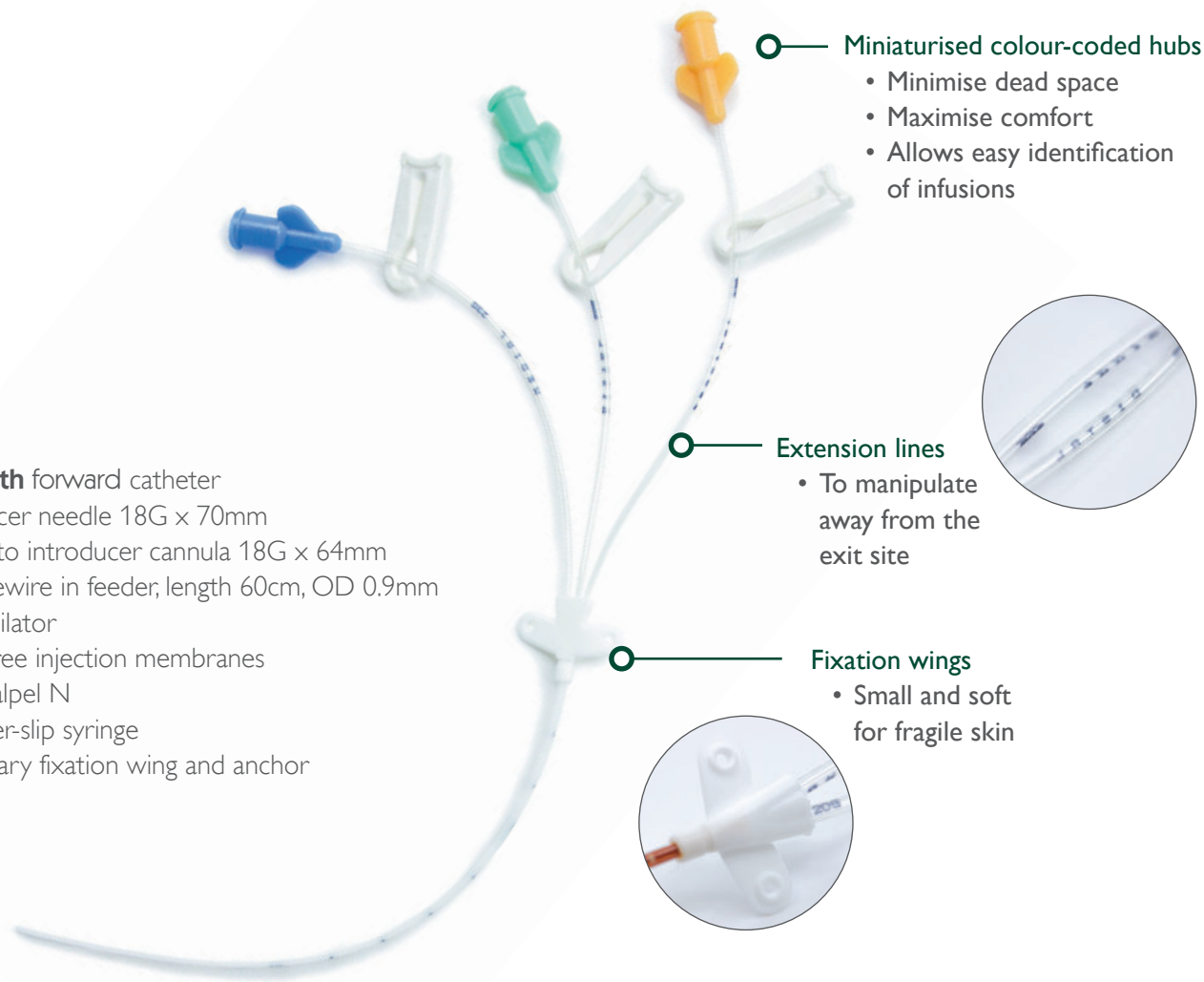

### Choice of catheter lengths and sizes:


Offering practitioners, a choice when selecting a catheter to suit the patient and therapy needs.

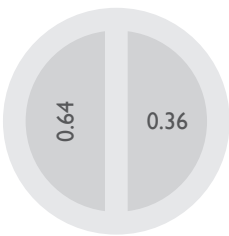


3Fr




4.5Fr

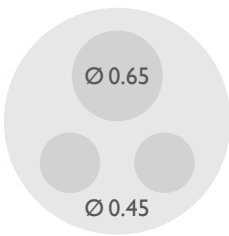

4.5Fr

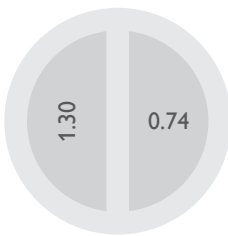

5.5Fr

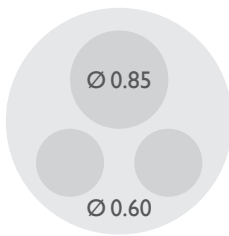

5.5Fr

multicath forward

Code	Lumen	Length (cm)	Catheter (PUR) Size (Fr)	OD (mm)	Flow rate & priming volume (mL/min)		
					Dist	Med	Prox
1202.063	2	6	4.5	1.5	22 0.20	-	24 0.20
1203.082	2	8	4.5	1.5	19 0.22	-	22 0.22
1203.123	2	12.5	4.5	1.5	15 0.25	-	18 0.25
157.508	2	8	5.5	1.85	46 0.34	-	77 0.38
1233.061	3	6	4.5	1.5	14.8 0.12	7.1 0.11	8.1 0.11
1233.081	3	8	4.5	1.5	14.4 0.14	6.7 0.12	7.5 0.12
1233.101	3	10	4.5	1.5	13.5 0.16	5.8 0.14	6.1 0.14
1231.121	3	12.5	4.5	1.5	10 0.18	3 0.16	4 0.16
155.167	3	16	7.5	2.5	63 0.38	34 0.33	36 0.38
155.207	3	20	7.5	2.7	60 0.45	30 0.34	30 0.36
0155.508	3	8cm	5.5fr	1.8	27ml /0.22	16ml / 0.21	17ml / 0.24
0155.512	3	12.5cm	5.5fr	1.8	24ml / 0.24	12ml/0.23	13ml / 0.23

smartmidline

smart**midline** is our CT rated midline (up to 325psi and 7mL/s). It is simple and quick to place using the direct Seldinger technique. Suitable for a wide range of clinical indications from 2Fr to 5Fr and 4cm to 25cm, it is also an integrated extension line, which helps to reduce the risk of post insertion complications.

Code	Size (Fr)	Length (cm)	Guidewire (cm)	Max flow rate* (mL/s)	Box qty
128.1208	2	8	26	1.0	20
128.1310	3	10	26	1.5	20
128.1315	3	15	40	1.0	20
128.1415	4	15	50	5.0	20

\*Viscosity 11.8cP



**CVC Webinar**  
 European Specialist in Veterinary Anaesthesia and Analgesia pinpoints potential mistakes and varying techniques of CVC placement in this video guide. Scan the QR code for more details.




# Nerve blocks

## echoplex needles

Designed for clinicians who use ultrasound and stimulation techniques to locate the peripheral nerve sheath.

Code	Size (G)	Length (mm)	Box qty
6194.103	21	100	20
6194.123	21	120	20
6194.153	21	150	20
6194.253	23	25	20
6194.353	23	35	20
6194.503	22	50	20
6194.853	21	85	20

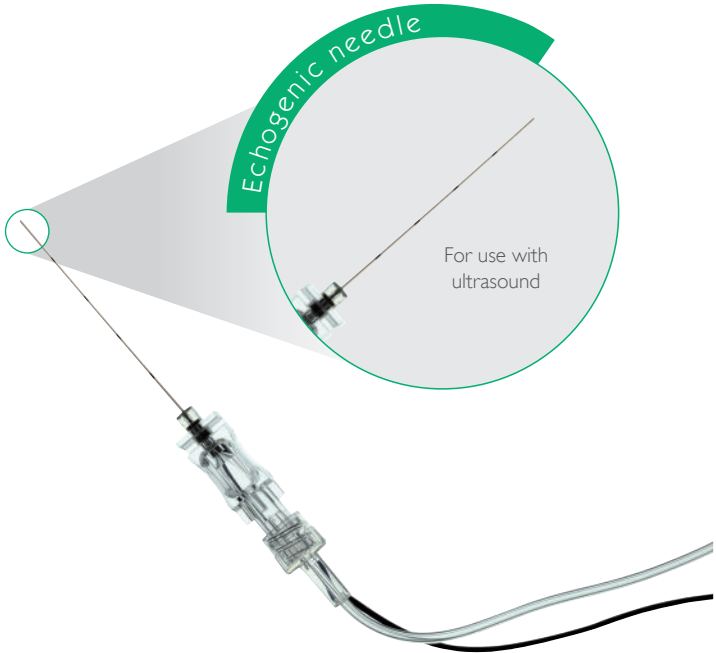

## plexygon

Nerve stimulator allowing nerve location through the emission of a low intensity electrical current.

Code	Description	Box qty
750131	plexygon nerve stimulator	1
08750302A	Replacement cable (Red and black)	1
3900.15	Aseptic handling pouch 12 x 40cm	80



“ We have experienced fewer problems with blood pressure, and we require lower levels of opiates. The dog can reliably go home the same day because the effects will last for around 8 to 12 hours. ”

**Mark Holmes**  
Clinical Director, Avondale Vets



Speak to a representative about a trialling a device. Various purchase options available.





Seeing a dog up and enjoying a meal not long after the operation rather than sitting salivating in their kennel was really positive and gave us a good indication immediately that nerve blocks were beneficial for patients.

It's an easy technique for the vets & nurses to learn, there really isn't too much to worry about with needle placement using the nerve stimulator. It was great for the nursing teams because they found they needed lower isoflurane levels during the surgery with the beneficial effect of keeping blood pressure levels more constant, with less fewer fluctuations in isoflurane levels needed. Everything just felt a bit calmer and smoother.

Introducing the Plexygon nerve stimulator with the training and support we received from Vygon made a lot of sense when we were able to see such nice clear-cut benefits for patients.



**Peter Haggis**

Advanced Practitioner in Small Animal Orthopaedics, Wilbury Vets



**Plexygon Nerve Stimulator Webinar**

European Specialist in Veterinary Anaesthesia and Analgesia pinpoints the benefits of using the Plexygon Nerve Stimulator in this video guide. Scan the QR code for more details.



# Cannula

## intraflon non-ported

A short IV cannula in PTFE with transparent ultra-thin walls to optimise flow rate.



Code	Colour	Size (G)	ID (mm)	OD (mm)	Length (cm)	Max flow rate* (mL/s)	Box qty
121.06	<div></div>	24	0.4	0.7	18	18	50
121.08	<div></div>	22	0.5	0.8	28	24	50
121.10		20	0.7	1.0	30	52	50
121.13	<div></div>	18	0.9	1.3	38	100	50
121.16	<div></div>	16	1.2	1.6	52	174	50
121.21	<div></div>	14	1.6	2.1	52	260	50
121.27	<div></div>	12	2.2	2.7	52	450	50

## intraflon non-ported longer

Code	Colour	Size (G)	ID (mm)	OD (mm)	Length (cm)	Max flow rate* (mL/s)	Box qty
122.10		20	0.7	1.0	36	50	50
122.13	<div></div>	18	0.9	1.3	47	95	50
122.16	<div></div>	16	1.2	1.6	60	170	50
122.21	<div></div>	14	1.6	2.1	80	250	50
122.27	<div></div>	12	2.2	2.7	80	430	50



## Haemodynamic Monitoring

Mostcare Up is haemodynamic monitoring at its most effective. With accurate patient monitoring minute by minute, beat to beat. With its intuitive, simple to use and customisable interface, you can set up and monitor a patient's cardiac output in 60 seconds, with no need for calibration.

Mostcare Up uses a robust patented algorithm known as the Pressure Recording Analytical Method (PRAM™), giving real time data to help you make the right decisions at the right time. Our endless system has unlimited patient usage and does not require consumables, ensuring you can optimise its usage without additional cost and helps improve sustainability. Mostcare Up can also connect to your patient information systems to help you stay more efficient.

Research article 'Briganti et al A preliminary study evaluating cardiac output measurement using Pressure Recording Analytical Method (PRAM) in anaesthetized dogs' available upon request.



### Fast

60 second setup



### Simple

No calibration required  
Intuitive, customisable interface



### Accurate

Robust PRAM™ algorithm  
Patented dynamic filter to enhance signal quality



### Lifetime value

Unlimited patient usage  
No consumables required, improving sustainability



### Data management

HL7 compatible to enable data transfer



Speak to a representative about a trialling a device.  
Various purchase options available.

## Procedure packs

A **vyset** procedure pack allows you to create a high quality pack to match your local protocol, incorporating all the items you need.

**vyset** procedure packs promote standardised insertion procedures, ensuring maximum barrier precautions and reduced risk of infection, ideally suited for many therapy areas including:

### The range

- Theatre & orthopaedic packs
- CVC packs
- Cannulation packs
- Biopsy packs
- Regional anaesthesia packs
- Soft tissue packs
- Arterial packs

“The custom packs are designed to our own specifications which is great for us. It saves us so much time all round and is no more expensive than buying everything separately. The boxes arrive sterilised and ready to use. We get them in as and when we need them and it's just a great service. As well as saving us time in ordering the products separately from different suppliers, we also save wear and tear on the autoclaves. We want to spend as much of our time with our patients and their owners as we can so the packs are just perfect for what we need.”

### Nicola Casement

Referral Veterinary Nurse, Manchester Veterinary Specialists







## Why choose vyset procedure packs?

From concept to implementation we will work with you to develop a high quality pack to meet your individual needs, giving you peace of mind throughout:

- Manufacturing
- Sterilisation
- Warehousing & distribution
- Implementation & training
- Stock management



Reduce waste



Save money



Reduce carbon footprint



Save on set up time

**Contact the vyset procedure pack team today to discuss your requirements**

t: 01793 746904

e: [veterinary@vygon.co.uk](mailto:veterinary@vygon.co.uk)


# Oncology

**MOSH**  
TESTED  
AND PASSED

## Qimono Qimooctopus

Qimooctopus with Qimo female and Qimo male is a multi-purpose catheter accessory, a double lumen IV extension tube with needleless connectors allowing intermittent injections or continuous infusions of fluids or medications.

Code	Box qty
7231.221	1

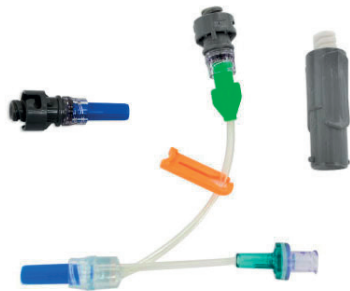

## Qimono Qimocap

Qimocap is a vial access cap fitted with Qimo female.

Code	Size (mm)	Box qty
7219.01	14	50
7219.02	20	25



## Qimono Qimospike V

Qimospike V is a vented spike for vials fitted with Qimo female.

Code	Box qty
7217.04	25



## Qimono Qimo male

Qimo male is a disinfectable closed needle-free safety connector.

Code	Box qty
7210.91	50



## Qimono Qimo female

Qimo female is a multi-purpose, closed, needle-free, I.V. connection device.

Code	Box qty
7210.02	50



## Non-sterile Chemotherapy gowns

Code	Description	Box qty
VHU69602	Non-sterile thumb loop Chemo gown L	100 (bagged in 10's)
VHU69603	Non-sterile thumb loop Chemo gown XL	100 (bagged in 10's)

## Sterile Chemotherapy gowns

Code	Description	Box qty
VHU44690	Aero chrome gown XL	28
VHU44693	Aero chrome gown L (Extra long)	30
VHU44694	Aero chrome gown XL (Extra long)	30

## Chemotherapy FFP3 facemask

Code	Description	Box qty
VHU62360	FFP3 Chemo mask with valve	20

## Chemotherapy gloves

Code	Description	Box qty
VHU50600M	Purple Nitrile-xtra glove XS	50
VHU50601M	Purple Nitrile-xtra glove S	50
VHU50602M	Purple Nitrile-xtra glove M	50
VHU50603M	Purple Nitrile-xtra glove L	50



# Enteral feeding



## nutrisafe2® PVC nasogastric feeding tubes

nutrisafe2 is a unique, non-luer system that is incompatible with other devices to prevent misconnections. It provides safe and secure nasogastric and naso-oesophageal delivery of vital nutrition. Includes colour-coded caps for easy identification of gauge size.

Code	Size (Fr)	Length (cm)	Box qty
DEHP-free PVC Feeding Tube – Short-term use			
0364042	4	40	25
0361042	4	50	50
0363042	4	75	25
0362042	4	125	25
0364052	5	40	25
0361052	5	50	50
0363052	5	75	25
0362052	5	125	25
0364062	6	40	25
0361062	6	50	50
0363062	6	75	25

Code	Size (Fr)	Length (cm)	Box qty
DEHP-free PVC Feeding Tube – Short-term use			
0362062	6	125	25
0364082	8	40	25
0361082	8	50	50
0363082	8	75	25
0362082	8	125	25
0361102	10	50	50
0363102	10	75	25
0362102	10	125	25
0363122	12	75	25
0362122	12	125	25

## nutrisafe2 PUR naso-oesophageal feeding tubes

The radiopaque feeding tubes come with size colour-coded cap and are suitable for short-term use. More sizes available on request.

Code	Size (Fr)	Length (cm)	Box qty
Polyurethane Feeding Tube – For use up to four weeks			
1362102	10	125	25
1363042	4	75	25

Code	Size (Fr)	Length (cm)	Box qty
Polyurethane Feeding Tube – For use up to four weeks			
1363082	8	75	25
1363122	12	75	25

## nutrisafe2 drawing-up straws

Drawing-up straws allow the filling of nutrisafe2 syringes with feed or medication from all sizes of containers and bottles.

Code	Length (cm)	ID (mm)	OD (mm)	Box qty
0817052	5	2.0	3.3	100
0817153	15	3.0	4.2	50



## nutrisafe2 syringes

nutrisafe2 syringes have a purple plunger to aid instant identification as enteral devices. For safety they only connect with the other items in the nutrisafe2 range. They are graduated and available in a wide range of sizes.

Code	Description	Box qty
1015014	nutrisafe2 1mL syringe	100
1015024	nutrisafe2 2.5mL syringe	100
1015054	nutrisafe2 5mL syringe	100
1015104	nutrisafe2 10mL syringe	100
1015214	nutrisafe2 20mL syringe	80
1015604	nutrisafe2 60mL syringe	60



## MIC-KEY G-tube kits\*

MIC-KEY G-tubes are designed for the delivery of enteral nutrition, medication and fluids to the stomach via a primary or pre-existing gastrostomy stoma.

Code	Size* (Fr)	Stoma length** (cm)	Balloon volume (mL)	Box qty
VKC8140-1225	12	2.5	3.5	1
VKC8140-1227	12	2.7	3.5	1
VKC8140-1230	12	3.0	3.5	1
VKC8140-1235	12	3.5	3.5	1
VKC8140-1240	12	4.0	3.5	1

\*Size range: 12Fr, 14Fr, 16Fr, 18Fr, 20Fr, 24Fr available  
 \*\*Stoma length range: 0.8cm to 5.0cm available



## MIC & MIC-KEY introducer kit\*\*\*

Redefining fluoroscopic, endoscopic and laparoscopic primary enteral feeding tube placement

Code	Description	Box qty
VKC98430-00	MIC & MIC-KEY G 12Fr Introducer Kit Includes 16Fr Dilator	1

\*\*\*Full MIC range also available. Please ask for details.



# Theatre

## Mini Red-O-Pack

120mL container in PE/PVC for drainage of small volumes after surgery.  
Supplied with Redon-type suction drain, Alene needle, and tape.

Code	Description	Suction drain size (Fr)	Container (mL)	Box qty
WKRAP0614101251	Mini Red-O-Pack	6	120	1
WKRAP0814101251	Mini Red-O-Pack	8	120	1



## Drape towels with adhesive

Two-layer general purpose sterile adhesive-edge drapes available in a variety of sizes.  
Adhesive edges provide an effective seal, preventing contamination of the surgical field from beneath the drape.

Code	Description	Size (cm)	Box qty
4126.03	Drape towels with adhesive	80x90	50
4126.05	Drape towels with adhesive	90x100	50
4129.05	Drape towels with adhesive	180x200	20



## Universal drape set

Set of four sterile adhesive-edged two-layer drapes suitable for many veterinary procedures including laparotomy. Ideal for small to medium sized veterinary patients.  
Suitable drape set for ventral midline ovariohysterectomy in the bitch.

Code	Description	Size (cm)	Box qty
4639.12	Adhesive-edged utility drapes x4	70x75	30

## Fenestrated drape towels

Code	Description	Size (cm)	Box qty
3445.45	Fenestrated towel with adhesive	45x45	80
3475.75	Fenestrated towel with adhesive	75x75	30
3475.90	Fenestrated towel with adhesive	75x90	30
3410.100	Fenestrated towel with adhesive	100x100	30

## Gowns

Sterile Surgical Gowns made from SMS material with hand towels, raglan sleeves, book folded for aseptic donning of gown.

Code	Description	Box qty
VHU11783	Sterile Surgical Gown Large with hand towels	40
VHU11784	Sterile Surgical Gown Extra Large with hand towels	40
VHU90002	Sterile Surgical Gown Small with hand towels	36
VHU90112	Sterile Reinforced surgical gown Large with hand towels	32
VHU90142	Sterile Reinforced Surgical gown extra large with hand towels	30
VHU90143	Sterile Reinforced Surgical gown Extra Large – Extra Long with hand towels	32





# Thoracic trocar and drain

Transparent radiopaque PVC thoracic trocar and drain with 2cm graduations and two lateral eyes. Handle ensures controlled insertion.

Code	Description	Size (Fr)	Length (cm)	Box qty
626.12	Thoracic trocar and drain	12	25	20
626.16	Thoracic trocar and drain	16	25	20
626.20	Thoracic trocar and drain	20	40	20
626.24	Thoracic trocar and drain	24	40	20
626.28	Thoracic trocar and drain	28	40	20





# Urinary

## Adjustable KatKath® catheter

All the benefits of the standard KatKath with an adjustable wing.

Code	Description	Box qty
IVTAAK3518K	Adjustable KatKath 3.5Fr/18cm	10



### KatKath benefits:

- **Simple to use** – unblock the cat. Leave in place
- **Non-reactive material** – low friction. Vital when dealing with an inflamed oedematous urethra
- **End hole** – allowing concurrent flushing during placement
- **Available with a rotating suture collar** – allowing natural rotation when in-situ
- **Kink and crush-resistant** – allowing placement of the catheter without the use of a stylet
- **Tapered distal tip** – allowing atraumatic catheter placement.

Code	Description	Box qty
1VTAAKK3511	Standard KatKath 3.5Fr/11cm	10
1VTAAKK3514	Standard KatKath 3.5Fr/14cm	10



# Aitecs infusion pumps

## aitecs 2017 syringe pump

The aitecs 2017 syringe pump is a configurable and versatile syringe pump offering enhanced safety in all clinical environments. The common user interface across the aitecs range promotes ease of use and simplifies device training needs. It can be used as a stand-alone device or be stacked together with other aitecs pumps on the aitecs infusion docking station providing connectivity compatible with the aitecs infusion management system, Patient Data Management Systems (PDMS) and Hospital Information Systems (HIS) via the open system architecture.

Code	Box qty
VMM2017TCI	1




## aitecs 3017 large volumetric pump

The aitecs 3017 large volumetric pump is a configurable multi-purpose infusion pump, with a common user interface to the aitecs range. It provides patient safety by facilitating use from basic infusions to advanced clinical applications. It can be used as a stand-alone device or be stacked together with other aitecs pumps on the aitecs infusion docking station providing connectivity compatible with the aitecs infusion management system, Patient Data Management Systems (PDMS) and Hospital Information Systems (HIS) via the open system architecture.

Code	Box qty
VMM3017TCI	1



We offer our pumps at no cost as the charges are offset against the consumables you purchase.

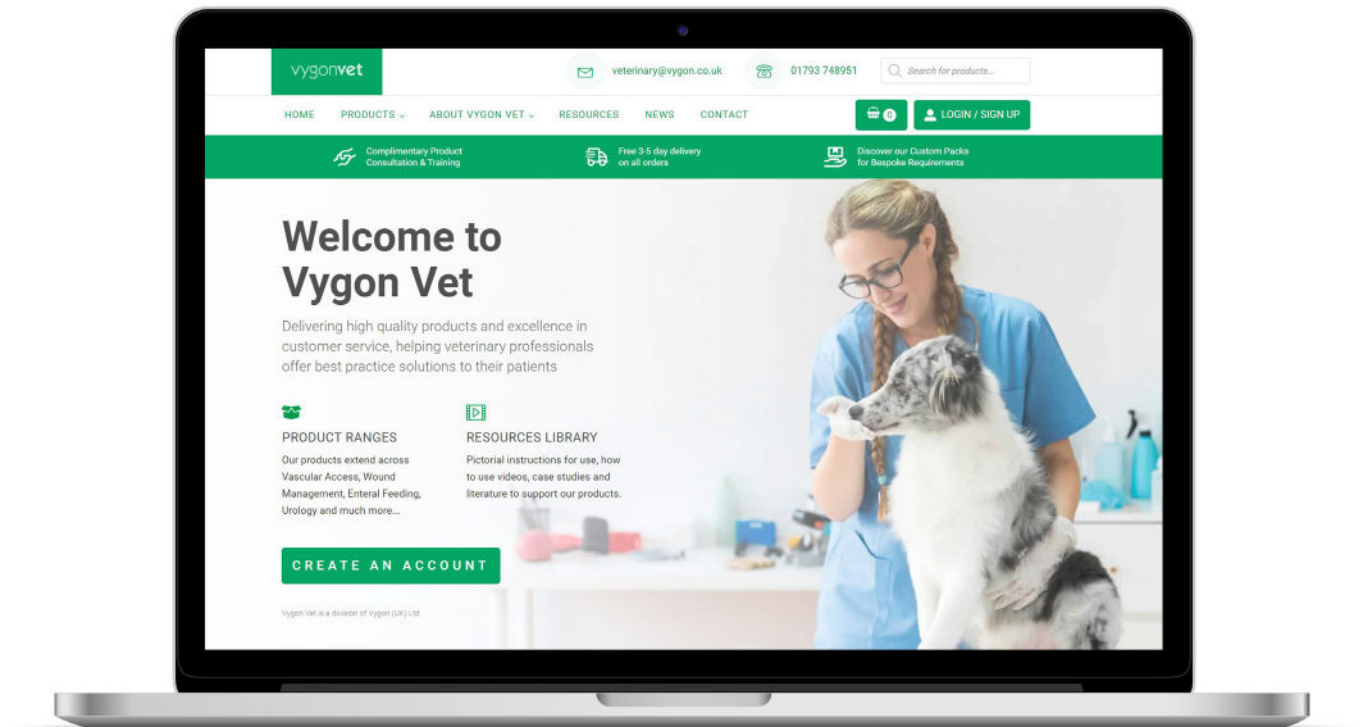

# Have you registered?

Sign up to an online account today and reap the benefits




## SAVE TIME

The site features all our products for instant online purchases.



## ACCOUNT SPECIFIC PRICING

Having an account witha Vygon Vet can save you up to 20% on certain products.

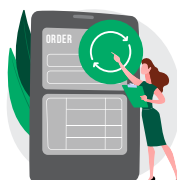

## REPEAT ORDERS WITH EASE

It takes just one click to repeat orders from previous purchases.

Take a look  [vygonvet.co.uk](https://vygonvet.co.uk)

## Our distributors

You can also buy direct from our distributors



[covetrus.co.uk](http://covetrus.co.uk)



[onlinestore.mwiah.co.uk](http://onlinestore.mwiah.co.uk)



[veterinarysales.co.uk](http://veterinarysales.co.uk)



[dmsveterinary.com](http://dmsveterinary.com)



[jakmarketing.co.uk](http://jakmarketing.co.uk)



[nvsweb.co.uk](http://nvsweb.co.uk)



[pioneervet.co.uk](http://pioneervet.co.uk)



[vet-direct.com](http://vet-direct.com)

FOR FURTHER INFORMATION, PLEASE CONTACT: [info@vygon.co.uk](mailto:info@vygon.co.uk)

The specifications shown in this leaflet are for information only and are not, under any circumstances, of a contractual nature. This brochure has been printed on responsibly sourced and sustainable material. To help us reduce our carbon footprint all of our literature is available electronically either from your Product Specialist, on our website or by emailing [info@vygon.co.uk](mailto:info@vygon.co.uk)

VYGON (UK) LTD, THE PIERRE SIMONET BUILDING, V PARK,  
GATEWAY NORTH, LATHAM ROAD, SWINDON, WILTSHIRE SN25 4DL  
RECEPTION: +44 (0)1793 748800 [WWW.VYGON.CO.UK](http://WWW.VYGON.CO.UK)

    [@vygonuk](https://www.instagram.com/vygonuk)

